

Product Description

Pioneering GTPase and Oncogene Product Development since 2010

CDC42(G12V) PROTEIN

产品名称: Cdc42(G12V) 突变蛋白

货号: 10109

产品全名: Cdc42 Protein G12V 突变蛋白

基因符号 Cell division cycle 42, G25K, CDC42Hs

Source: Human, recombinant full length, His6-tag

Expression 种属反应性: E. col

分子量: 21 kDa

纯化: >95% by SDS-PAGE

Introduction: Small GTPases are a super-family of cellular signaling regulators. Cdc42 belongs to the Rho sub-family of GTPases that regulate cell motility, cell division, and gene transcription. GTP binding increases the activity of Cdc42, and the hydrolysis of GTP to GDP renders it inactive. GTP hydrolysis is aided by GTPase activating proteins (GAPs), while exchange of GDP for GTP is facilitated by guanine nucleotide exchange factors (GEFs).

Amino Acid Sequence (1-191, G12V)

MQTIKCVVVGDVAVGKTCLLISYTTNKFPSEYVPTVFDNYAVTVMIGGEPYTLGLFDTAGQEDYDRL RPLSYPQTDVFLVCFSVVSPSSFENVKEKWVPEITHHCPKTPFLLVGTQIDLRDDPSTIEKLAKNKQ KPITPETAEKLARDLKAVKYVECSALTQKGLKNVFDEAILAALEPPEPKKSRRCVLL

Properties

Physical Appearance (form): Dissolved in 20mM Tris-HCl, pH8.0, 150mM NaCl.


Physical Appearance (form): White or clear

Concentration: 1 mg/mL

Storage: -80°C

Preparation Instructions:

Centrifuge the vial before open the cap and reconstitute in water. Adding of 10 mM β -mercaptoethanol or 1 mM DTT into the solution to protect the protein is recommended and using of non-ionic detergents such as n-Dodecyl β -D-maltoside (DoDM) or polyethylene detergents (e.g. C12E10) also help to stabilize the protein. Avoid repeated freezing and thawing after reconstitution. The purity of His-tagged Cdc42 G12V was determined by SDS-PAGE and Coomassie Brilliant Blue Staining.

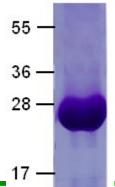



Product Description

Pioneering GTPase and Oncogene Product Development since 2010

References:

- 1. Garrett, W. S. . et al., Cell 102: 325-334, 2000.
- 2. Irie, F. et al., Nature Neurosci. 5: 1117-1118, 2002.
- 3. Kawasaki, Y. et al., Oncogene 26: 7620-7627, 2007.
- 4. Manser, E. et al., Nature 363: 364-367, 1993.
- 5. Musch, A. et al., EMBO J. 20: 2171-2179, 2001.
- 6. Nalbant, P. et al., Science 305: 1615-1619, 2004.
- 7. Shen, Y. et al., Dev. Cell 14: 342-353, 2008.
- 8. Wu, W. J. et al., Nature 405: 800-804, 2000.
- 9. Wu, X. et al., Genes Dev. 20: 571-585, 2006.
- 10. Zheng, Y. et al., J. Biol. Chem. 271: 33169-33172, 1996.